

Milli-Q



Milli-Q® HX 7040-7150 / HR 7060-7220

User Manual



Contents

Legal Information	5
Safety Information	6
Recycling	7
Preface	8
Specifications	9
Feed Water Specifications.....	10
Product Water Specifications.....	10
System Dimensions.....	12
Weight Specifications.....	13
Electrical Specifications.....	13
Environmental Specifications.....	14
System Overview	15
Operating Principle	16
Distribution loop management and configurations.....	17
System Display	18
How to Navigate Between Overview and Workspace Screens.....	18
Overview Screen.....	19
How to Change the System from Standby to Ready.....	23
Workspace Screens.....	25
Glance Workspace.....	26
Maintenance Workspace.....	30
Configuration Workspace.....	32
How to View Hidden Applications.....	35
Alarms and Alerts - Definition.....	36
How to Acknowledge Alarms and Alerts.....	37
Screensaver.....	38
Maintenance	39
How to Use Maintenance Wizard.....	40
RO Cleaning tools.....	41
RO Cleaning.....	42
Communication Ports & Software	43
How to Change the Network Configuration.....	44
How to Power On and Power Off the Water System	46
How to Use the Distribution Schedule	47
How to Flush the Tank	49

Ordering Information..... 50
 Consumable Catalogue Numbers..... 51
 Accessory Catalogue Numbers 52
 System Catalogue Numbers..... 53

Appendix..... 54
 Display Icon Description..... 55
 System Operating modes..... 58
 Emergency Procedure in Multi-system configuration..... 59

Legal Information

Notice

It has always been Millipore S.A.S. policy to continuously improve its products.

The information in this document is subject to change without notice and should not be construed as a commitment by Millipore S.A.S. Millipore S.A.S. assumes no responsibility for any errors that might appear in this document. This manual is believed to be complete and accurate at the time of publication. In no event shall Millipore S.A.S. be liable for incidental or consequential damages in connection with or arising from the use of this manual.

We manufacture and sell water purification systems designed to produce pure or ultrapure water with specific characteristics ($\mu\text{S/cm}$, $^{\circ}\text{C}$, TOC, CFU/ml, EU/ml) when it leaves the water purification system provided that the systems are fed with water quality within specifications, and properly maintained as recommended by the supplier.

We do not warrant these systems for any specific application. It is up to the end user to determine if the quality of the water produced by our systems matches his expectations, fits with norms/legal requirements and to bear responsibility resulting from the usage of the water.

Product Warranty and limitation of liability

The applicable warranty and limitation of liability for the products listed in this publication may be found at <http://www.merckmillipore.com/> or <http://www.emdmillipore.com/> within the "Conditions of Sale" applicable to your purchase transaction.

Copyright

© 2017 MILLIPORE S.A.S. PRINTED IN FRANCE. ALL RIGHTS RESERVED. THIS BOOK OR PARTS THEREOF MAY NOT BE REPRODUCED IN ANY FORM WITHOUT THE WRITTEN PERMISSION OF MILLIPORE S.A.S.

The photographs illustrating the products are non-contractual.

Trademarks

Elix and Progard are registered trademarks of Merck KGaA, Darmstadt, Germany.






The " M " logo is registered trademark of Merck KGaA, Darmstadt, Germany.


All other trademarks are trademarks of their respective manufacturers.

Document Reference: USER-LRG2-ELIX-EN

Revision: V4.0

Safety Information

Symbol	What it means
	This <u>UV RADIATION</u> sticker is used to refer to a position on the system cabinet or inside of it where exposure to UV light is possible.
	This <u>HAZARD</u> sticker is used to refer to a position on the system cabinet or inside of it that could be hazardous.
	This <u>ELECTRICAL GROUND</u> sticker is used to refer to a position on the system cabinet or inside where an electrical ground connection is made.
	This <u>ELECTRICAL HAZARD</u> sticker is used to refer to a position on the system cabinet or inside where an electrical danger could exist.
	This <u>HOT SURFACE</u> sticker is used to refer to a position on the system cabinet or inside where a hot surface could exist.

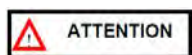

- The System must be connected to a source of electrical power that is earth grounded.
- Only an authorized person following the local safety regulations can work on this equipment.
- Unplug the electrical power cord before plugging in or removing any components on the electronic board(s).

Recycling

Directive 2002/96 EC: For European users only.



The symbol "crossed bin" on a product or its packaging indicates that the product should not be treated like household waste when discarded. Instead, the product should be disposed of at a location that handles discarded electric or electronic equipment.

Proper disposal of equipment containing electric or electronic components will help to reduce pollution effects to the environment or to human health. Proper recycling of these products helps in environmental preservation and helps to protect natural resources. For more information about recycling of products containing electric or electronic components, please contact your local recycling representative or organization.

Preface

Thank you for purchasing our water purification system.

For correct operation, read and fully understand the contents of this user manual before attempting to use the system. It is recommended to store this user manual in a safe and convenient place where it can be easily referred to when required.

This user manual is intended for use with a Milli-Q[®] HX 7040-7150 - HR 7060-7220 water purification system.

The word "system" is used to refer to Milli-Q[®] HX 7040-7150 - HR 7060-7220 unless specified otherwise in the text.

The abbreviations "LC" and "HC" used after HX 7040-7080 and HR 7060-7120 refer to systems for Low Chlorine and High Chlorine feed water.

The words "display" or "MMI" (Man-Machine Interface) are used to refer to the user interface in this document.

The words "master" and "slave" have been replaced respectively by "primary" and "auxiliary" for systems with software package >35.

Depending on the system type and options that are purchased, some of the features described may not apply to the system you are using.

For any questions or requests, please contact your Sales Specialist or Qualified Service Representative.

Specifications

Feed Water Specifications

The system has been designed to operate within feed water specifications:

Parameter	Value or range (HX and HR Systems)
Pressure	2 – 6 bar
Flow rate	> 10 L/min @ 2 bar
Feed water type	Potable water
Temperature	5 – 35° C
Conductivity	10 – 2000 µS/cm @ 25° C
pH	4 - 10
Hardness (as CaCO ₃)	< 300 ppm
Silica concentration	< 30 ppm
Carbon Dioxide concentration (CO ₂)	< 30 ppm
Langelier Saturation Index (LSI)	< 0.3
Fouling Index (FI ₅) or Silt Density Index (SDI)	7 (*)
Total Organic Carbon (TOC)	1 ppm
Free chlorine for LC systems	< 1.5 ppm
Free chlorine for HC systems	1.5 ppm - 3 ppm

* < 12 when the optional UF pretreatment unit is installed.

Product Water Specifications

The system has been designed to produce water according to specifications when operating within Feed Water Specifications.

Parameter	HX Systems: Value or range	HR Systems: Value or range
Resistivity	> 5 MΩ.cm @ 25°C	NA
Conductivity	< 0.2 µS/cm @ 25 °C	95% ionic rejection (99% Particulates rejection)
Total Organic Carbon (TOC)	< 30 ppb	99% organic rejection for MW>200 Dalton
Bacteria	< 10 CFU/mL	NA
Silica	<3ppb (rejection > 99,9%)	NA

Applicable standards TO HX system:




The system has been designed to be compliant with water quality standards in the laboratory.

The system produces water that meets the following requirements:

- ISO 3696: 1996 Grade 2 Water for Analytical laboratory use
- GB/T6682-2008 Grade 2 Water for Analytical laboratory use
- ASTM D1193 type 2 (2006 Reapproved 2011) Reagent Grade Water
- Japanese Industrial Standard JIS K 0557, A3 – 2008
- European Pharmacopeia Purified water 8.0
- United States Pharmacopeia Purified water (USP37)
- Chinese Pharmacopeia (2010 appendix XVII A-227) Water for Pharmaceutical Purposes
- Japanese Pharmacopeia (17-2016) Purified Water

System Dimensions

HX and HR System dimensions:

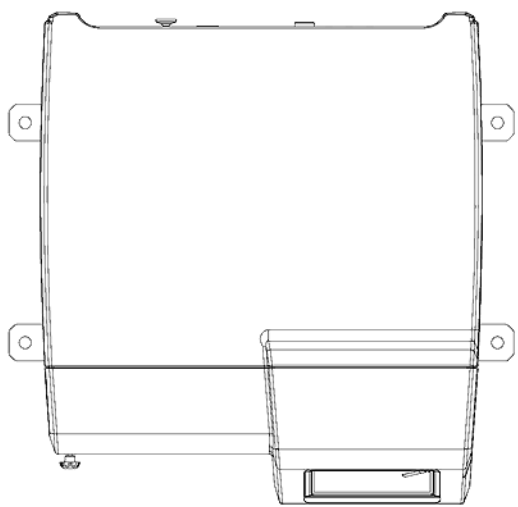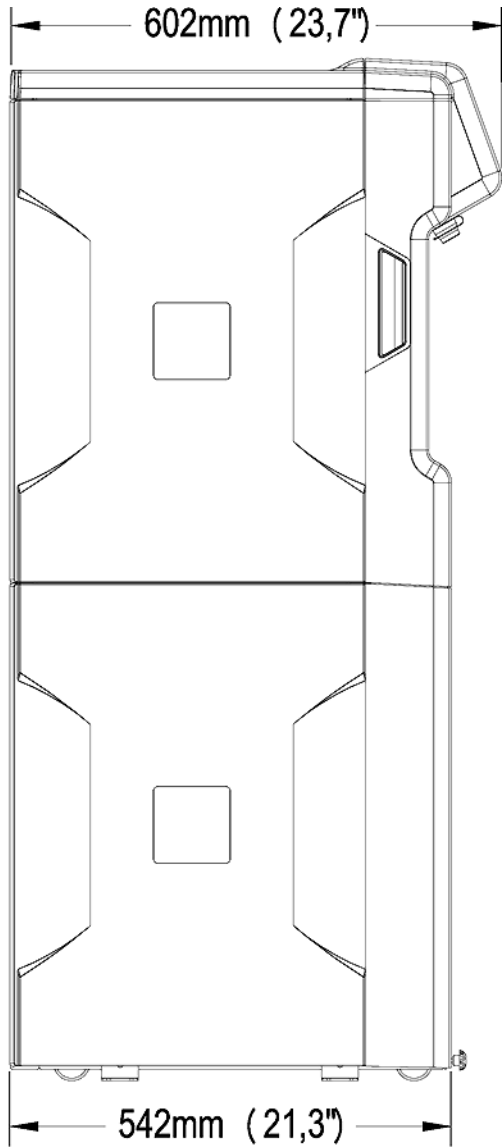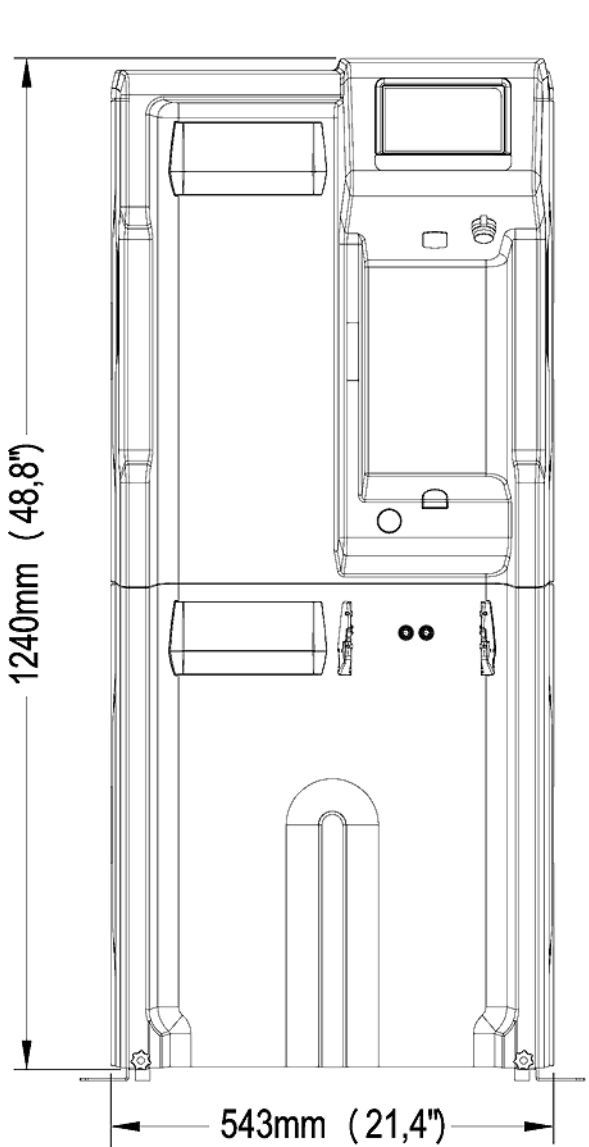

Weight Specifications

The location where the Water Purification System is installed will need to fully support its Operating Weight:

System type	Dry weight (Kg/Lb)	Shipping weight (Kg/Lb)
Milli-Q [®] HX 7040	78 / 172	97 / 213
Milli-Q [®] HX 7080	86 / 190	105 / 231
Milli-Q [®] HX 7120	94 / 208	113 / 249
Milli-Q [®] HX 7150	105 / 232	124 / 273
Milli-Q [®] HR 7060	72 / 166	91 / 207
Milli-Q [®] HR 7120	75 / 179	94 / 220
Milli-Q [®] HR 7170	78 / 192	97 / 233
Milli-Q [®] HR 7220	84 / 221	103 / 252

Shipping Weight is defined as a dry system in its shipping container. It does not include the consumable or any accessories.

Dry Weight is defined as a system without its shipping container. It does not include the consumable or any accessories.

The UF pre filtration option weight is 22 kg / 121.3 lb

Electrical Specifications

The systems can be supplied by Mains electricity 90-253 VAC with a frequency range from 48 to 62Hz.

System type	Voltage	Power Consumption (VA)
Milli-Q [®] HX 7040/7080	220-240 VAC @ 50/60 Hz	750
Milli-Q [®] HX 7120/7150	220-240 VAC @ 50/60 Hz	870
Milli-Q [®] HX 7040/7080	120 VAC @ 60 Hz	775
Milli-Q [®] HX 7120/7150	120 VAC @ 60 Hz	900
Milli-Q [®] HX 7040/7080	100 VAC @ 50/60 Hz	775
Milli-Q [®] HX 7120/7150	100 VAC @ 50/60 Hz	900
Milli-Q [®] HR 7060/7120	220-240 VAC @ 50/60 Hz	620
Milli-Q [®] HR 7170/7220	220-240 VAC @ 50/60 Hz	750
Milli-Q [®] HR 7060/7120	120 VAC @ 60 Hz	620
Milli-Q [®] HR 7170/7220	120 VAC @ 60 Hz	750
Milli-Q [®] HR 7060/7120	100 VAC @ 50/60 Hz	620
Milli-Q [®] HR7170/7220	100 VAC @ 50/60 Hz	750

The source of electrical power must be earth grounded.

Environmental Specifications


Environmental specifications have been defined for normal system operation.

Altitude	2000 m
Ambient operating temperature	10 – 40 °C
Relative humidity	80% up to 31°C (decreasing linearly to 50 % relative humidity at 40 °C)
Ambient storage temperature	0 – 40 °C
Ambient storage humidity	10 – 95 %
Installation category	I
Pollution degree	2
Noise level	< 50 db at 1 meter

For indoor use only

System Overview

System, components and connections.

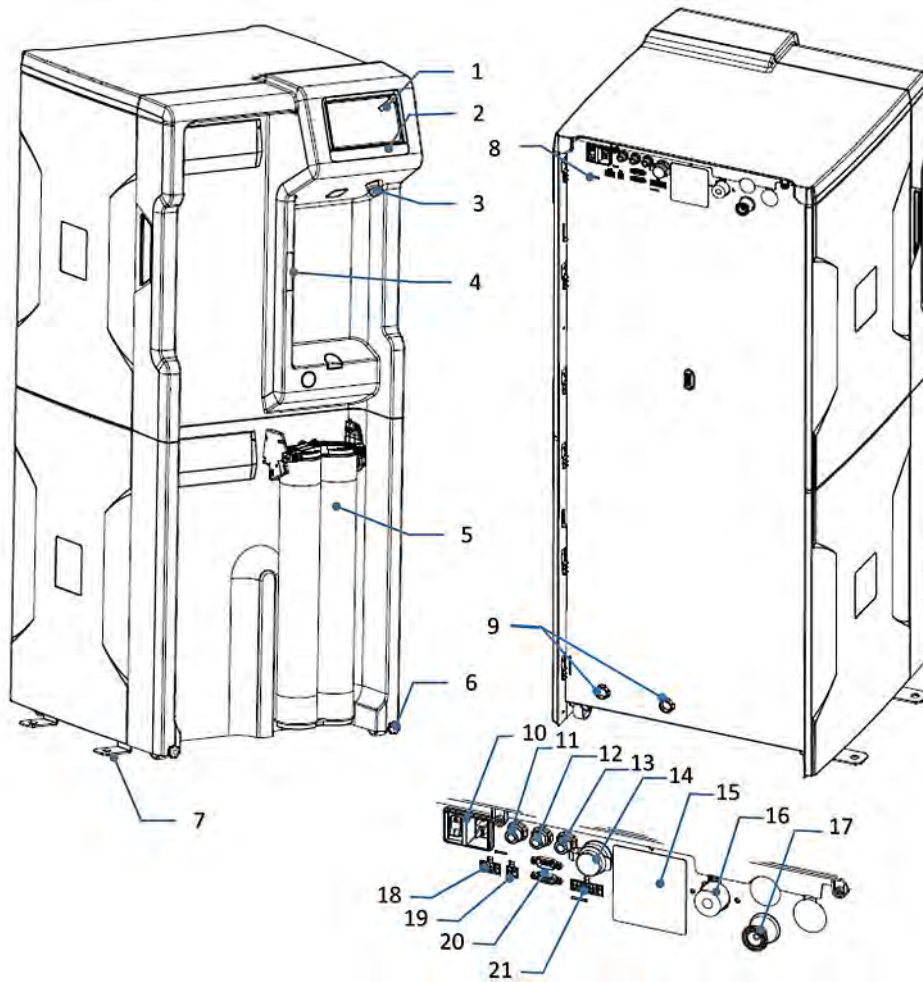

1. Display	12. DRP port (Data report: 2 Alarm output and 2 measure output)
2. System name	13. External solenoid valve port
3. USB port	14. Ethernet port
4. Matrix code (containing system serial number and type)	15. System type identification label
5. Progard [®] pack(s)	16. Inlet connector (3/4" BSP)
6. Wheel locking nut	17. Outlet connector (3/4" Tri-Clamp)
7. Securing bracket	18. External pretreatment port
8. Water system backpanel	19. Leak detector port
9. Waste outlet	20. Tank monitoring port (Tank level, overflow)
10. Power Entry (Mains) & switch	21. Device control port (ASM, UV lamp)
11. C2 Out port (Distribution interfaces)	




Note: The system covers can be cleaned or disinfected by using alcohol (Ethanol and Isopropanol).

Operating Principle


The Elix water system is designed to produce Type 2 water to feed a SDS 500 reservoir or a custom tank and to manage a distribution loop and accessories when installed. The system can be installed in a single or multi-system configuration.

- The water system purifies tap water using a Progard[®] Pack(s) and Reverse Osmosis followed by Elix[®] Technology. This process is called **Makeup (I)**.
- The purified water is then stored in a tank (SDS 500 or custom tank). This is referred to as **Storage (II)**.
- The stored water can be distributed to a loop. This process is called **Distribution (III)**.

HX 7040 flow schematic:

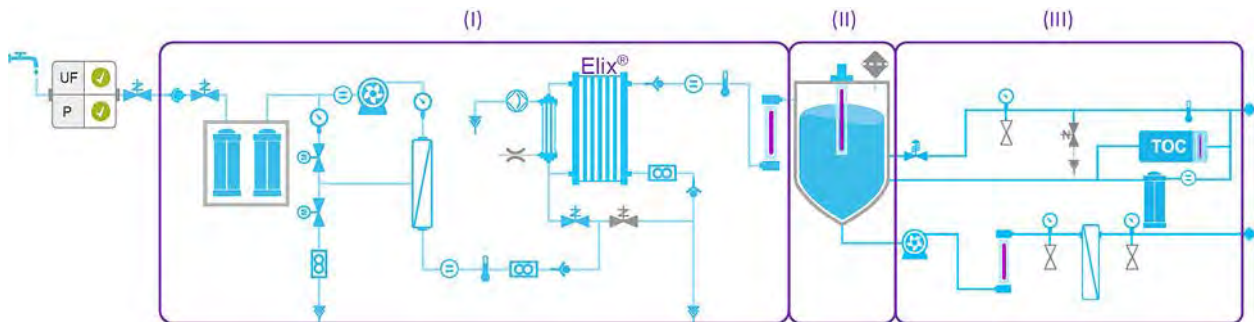

HR 7060 flow schematic:

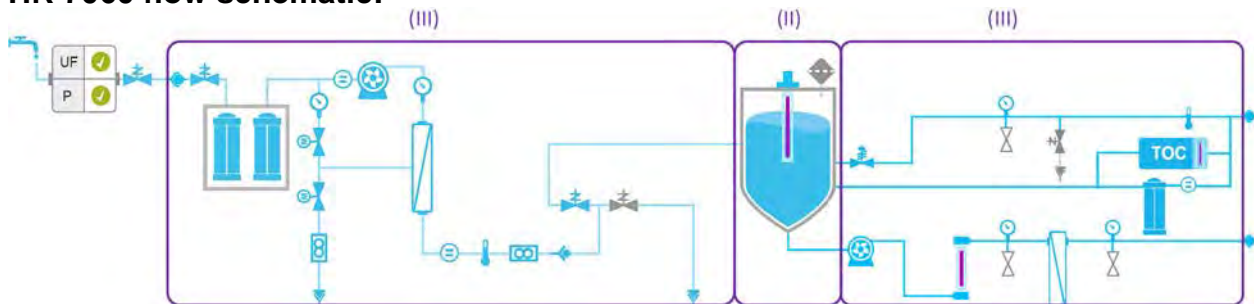

This view and more information about components can be retrieved using the **Flow**

Schematic application  at   
(Glance Workspace > Flow Schematic)

The volume of water in the Tank controls Makeup water production.

The storage part of the system is passive. The volume of water in the tank is a function of both the Makeup and Distribution processes.

Ready mode is used to get the processes operating.

Standby mode is used to stop the corresponding process in the event of a leak or maintenance purpose.

System operating modes are listed in the Appendix.

The Makeup and Distribution are two independent processes. When one of the processes is in standby, alarm stop or maintenance it will not stop the other process from working.

Distribution loop management and configurations

A system can be installed as a single system or as part of a multi-system configuration. If installed, a system can manage a distribution loop and its accessories.

System not managing a distribution loop

When the system does not manage a distribution loop, the system needs to get information from the external sensors as the system needs to stop in certain conditions, for example if an external tank is full or in the case of water detected. The following signals are taken into account by the system:

- Tank level
- Water Sensor
- Overflow

System in a single-system configuration managing a distribution loop

The system is connected to the external sensors and actuators linked to the distribution loop:

- The pump in simplex or duplex mode
- Tank level management
- The ASM
- UV lamp
- Dumping valve
- Leak detection
- Overflow

System in a multi-system configuration

In a multi-system configuration, the systems are connected together in order to share the same information (For example: tank level sensor information, or a water detection alarm that will stop all systems) and to synchronize the production processes.

In this configuration, one system is called the “Primary” and the other system(s) is/are called the “auxiliary(s)”. A maximum of 3 systems can be installed in this way.

The connections related to the distribution are connected to the primary system. In this configuration the primary system manages the distribution loop and the related consumable lifetimes.

System Display

The main display on the system is a touchscreen. Tapping display button icons moves the display between screens and launches applications.



The Overview screen is the main screen. The Workspace screens (Glance, Maintenance and Configuration) contain the system applications.

How to Navigate Between Overview and Workspace Screens

Procedure

- Navigation between the **Overview** screen and the **Workspace** screens is done using the display icon buttons on the bottom right of the screens.




- From **Overview** when tapping the arrow icon , **Glance Workspace** is always the first screen opened.
- The return to **Overview** can be made from any of the three **Workspace** screens by using the Exit icon .

Overview Screen

The **Overview** screen is the default view on the system display. It is divided into 3 sections and each section represents one of the 3 parts of the water system.

Figure 1: Overview Screen

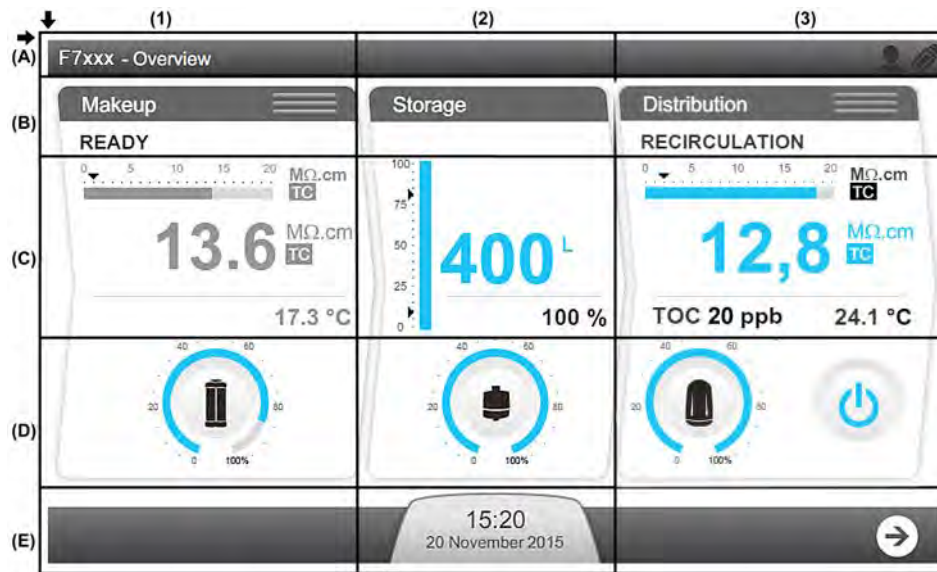

(1): Makeup - information on Makeup process

(2): Storage - information on Storage portion

(3): Distribution - information on Distribution process

(A): Top banner - System Information

A1: Screen Identifier - Overview

A3: Connectivity (external interfaces: Ethernet and USB)

(B): Water System Process Identifier and Operating Mode

B1: Makeup process is in **READY** mode.

B3: Distribution process is in **RECIRCULATION** mode. The Distribution process can also be **NOT ACTIVE** or in **MULTI-SYSTEM auxiliary** mode.


The system operating modes are described in the appendix.

(C): Measure Widgets (measure widgets are blue when the process is active or grey if it is not active)

C1: Measures from the Makeup process








	1	Bar graph showing water quality in MΩ.cm @ 25°C and the alarm setpoint (black triangle). Available only on HX systems.
	2	Displays Makeup water quality in a user defined unit. TC is displayed when Temperature Compensation (TC) mode is on. On HR systems, the unit is μS. If the water quality measured is under the setpoint, this value blinks between blue and red.
	3	Displays Makeup water temperature in a user defined unit.

C2: Measures from the Storage process



	1	<p>Bar graph showing tank level in %.</p> <ul style="list-style-type: none"> • Top arrow indicates Tank Filling Restart setpoint • Bottom arrow indicates TANK EMPTY alarm setpoint
	2	Displays volume of water in the Tank in a user defined unit.
	3	Displays tank level in %.

C3: Distribution process.

The distribution can be in the following modes:


	<p>Distribution in Forced mode:</p> <p>Distribution process is forced to ON. Distribution loop in this mode is always recirculating.</p> <p> Note: This function can be used to refresh the water from the distribution loop, for example after a loop sanitization or prior a sampling procedure use.</p>
	<p>Distribution in Scheduled mode:</p> <p>The Distribution is activated based on defined schedule.</p> <p> Note: It is recommended to activate this mode for normal use of the system and distribution loop.</p>
	<p>Distribution in Auto-recirculation mode:</p> <p>An automatic recirculation occurs once every hour (15 to 60 min/h depending on the software configuration).</p> <p> Note: It is recommended to activate this mode when the System and distribution loop are not used for a period of time.</p>
	<p>Distribution is Inactive :</p> <p>The system is not managing a distribution loop or the distribution loop is managed by a primary.</p>

The distribution panel will be displayed according to the distribution configuration of the system:

With distribution kits installed	Without distribution kits installed
	

With distribution kits installed:


To access the distribution modes, click on distribution mode button.

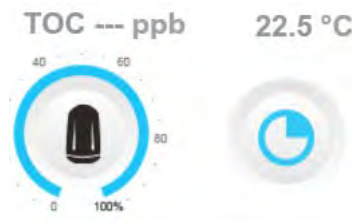



The distribution panels will show the three distribution modes. Select the desired mode:

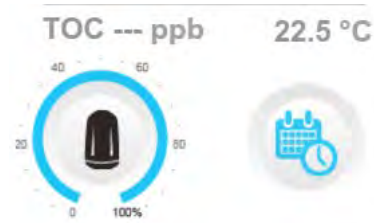

The distribution panel will return to the initial distribution screen:



(D): Consumable Status

Capacity remaining (%)	Consumable due for replacement (Blinking)
	


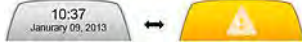

D1: Progard® status and gauge (Depending on system type, one or two Progard® can be displayed)

D2: Vent filter status and gauge

D3: Loop filter status and gauge

(E): Notification and Navigation bar

E2: Notification button showing date and time. Colour changes depending on system status.

No Alert(s) or Alarm(s) Present	Alert(s) Present (Blinking)	Alarm(s) Present (Blinking)
		

E3: Navigation to Workspace screen button.


How to Change the System from Standby to Ready

Keeping the system processes in Ready means the system automatically changes software modes to produce water when required.

Procedure

- !** **Important:** Standby should only be used to stop the corresponding process in the event of a leak or maintenance purpose.

If a process is not in a **READY** mode the operating mode will display **STANDBY**.

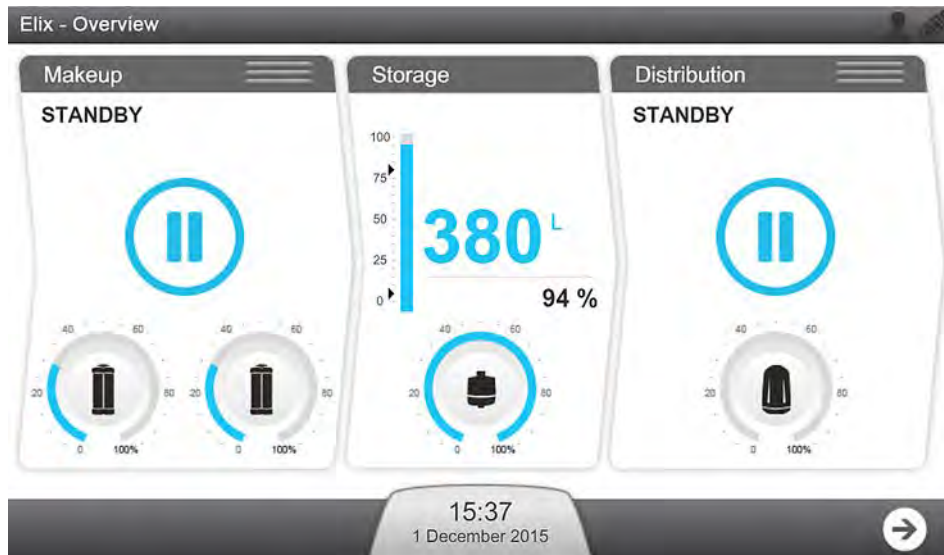

1. The Makeup and Distribution processes can be changed between **STANDBY** and **READY** using the three horizontal bars at the top of their respective window.


(I): Makeup; (II): Distribution



Note: When the system is not managing a distribution loop, the distribution panel is not active and in this case does not display the three horizontal bars.

2. Press on the icon on the horizontal bar of the process that you wish to change.

By clicking on the bar => Process mode switches:

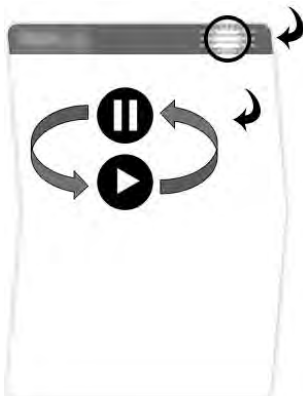

3. Confirm this action.



Once confirmed, the corresponding process enters **"PLEASE WAIT"**



Results

When both Makeup and Distribution processes are in Ready, the filling of the tank and loop management is automatic.

For more information, please refer to System Operating Modes in the Appendix.


Workspace Screens

There are three Workspace screens, **Glance**, **Maintenance** and **Configuration** which list all available applications.

Glance Workspace

The **Glance Workspace** contains applications providing system information.



 **Note:** When switching from an application back to the workspace screen, the last workspace screen selected will be opened.



Information Application gives information about your system.

Product Information: Unique manufacturing identifier of your system. You may need to provide some or all of this information when contacting your Sales or Service representative.

Software versions: Your water system is made of several electronic boards using different software. You may need to provide some or all of this information to your Service representative in case of troubleshooting. Software and firmware board versions are listed here.

Customer information: Details such as your location and address information can be seen here.



Service Application gives information related to Service.

Key contacts: Information about support key contacts is seen here (technical service representative or sales support contact information).

Service information: Information about your system's service agreements. It contains the installation details, and if applicable the contract name and the operational or preventive visit dates.



Consumables Application provides information related to your system's consumables.

Progard® pack: Lists the Progard® pack(s) details. The catalog number, installation date, processed volume and replacement date.


RO cartridge: Identifies your RO cartridge using type, lot and catalog number. The installation date can also be checked here.

Vent filter: Lists the tank vent filter type, lot and catalog numbers with the installation and replacement dates.

Loop filter: Lists your Loop filter details. The catalog number, installation date, processed volume and replacement date.



Measures Application gives information on the measurements from your system.

Water quality measurements:  provides information about each stage of the water purification process in the system.

- **Tap water:** Conductivity and temperature are monitored.
- **RO stage:** Conductivity, temperature and the permeate and ionic rejection outputs are monitored.
- **ELIX® stage:** On HX systems, the Elix® product conductivity and temperature are monitored.
- **Distribution stage:** If the kits are installed, TOC, Distribution Resistivity and temperature are displayed.

Hydraulic measurements & actuator statuses:  shows all hydraulic measurements and gives the status of each actuator.




- **Tap water:** The tap feed pressure and if applicable, the external solenoid valve status is displayed.
- **RO stage:** All parameters related to the RO stage are displayed. For example, the RO pressure, flowrate and recirculation and the different solenoid valve statuses.
- **Elix® stage:** On HX systems, the parameters related to the Elix® stage are displayed. In this application, the Elix® status, product flowrate and recovery, and the UV lamp status are monitored. TOC, distribution resistivity and temperature are also monitored if the options are installed in the system.
- **Storage stage:** The tank level, and, if applicable, the ASM UV lamp status are displayed.
- **Distribution stage:** If applicable, the distribution pump status and the dumping valve status are displayed.

Electrical measurements:  shows electrical measurements.

Electrical values related to the different system stages are monitored within the application.



Flow Schematic Application shows an overview of your system with a realtime updated flow schematic.

All system actuators are shown on the flow schematic. All of the different measurements are displayed, the values displayed can be sorted by water quality , hydraulic  or electrical measurement .

Pressing on a component or a consumable item will bring up a **tool tip**. The tool tip gives the item name and related information.



Components Application provides information about the major components of the system.

Pumps: Reverse Osmosis, Distribution and if applicable, Degasser pump details. The catalog number, installation and replacement date are displayed.

UV lamps: Catalog number, installation and replacement date are displayed for the following components: Makeup on HX systems only, and if installed, Automatic Sanitization Module, Distribution UV Lamp and TOC UV lamp.



Sanitization and Cleaning Application shows information about the sanitization and cleaning that is performed on the system. When the due date of a cleaning or sanitization nears, the system will display an alert.



Note: You can only visualize the consumable status information from this application. To perform consumable replacement, use the **Maintenance Workspace > Consumables Application**.

RO membrane CL₂ cleaning: The last and next CL₂ cleaning to be performed are displayed. If the option has been activated in the system.

RO membrane pH cleaning: The last and next pH cleaning to be performed are displayed. If the option has been activated in the system.



History Application provides historical information relating to your system. The System data can be filtered by date and the selected data can be exported in xml files.

Daily measurements: Historical daily system measurements and internal parameters can be viewed and retrieved. These parameters represent an indicator of the system performance.

Daily operations: The system creates a daily journal that contains information about all operations performed on it. Typical values registered include the water volume processed and the number of working hours.

Event record: Event records like power on, power off, operating modes or the operations that trigger alarms or alerts can be viewed and retrieved.



User Manual Application allows the system user manual to be downloaded (using USB or Ethernet connection).



Distribution Schedule Application shows the Distribution schedule (if the system manages a distribution loop).

Auto-recirculation: distribution loop automatically recirculates once per hour during the period scheduled. (Auto-recirculation cycle duration can be adjusted from 15 to 60min / hour by your Qualified Service representative)

Recirculation: distribution loop recirculation is forced during the period scheduled.

Tank flush: water is discarded from the reservoir during the period scheduled.



Diagnostic Application allows the system log files to be downloaded in order to be sent to Service Representative for diagnostic purposes.



Distribution Schedule Application shows the Distribution schedule (if the system manages a distribution loop).

Auto-recirculation: distribution loop automatically recirculates once per hour during the period scheduled. (Auto-recirculation cycle duration can be adjusted from 15 to 60min / hour by your Qualified Service representative)

Recirculation: distribution loop recirculation is forced during the period scheduled.

Tank flush: water is discarded from the reservoir during the period scheduled.

Maintenance Workspace

The **Maintenance Workspace** contains applications enabling maintenance and cleanings to be performed.

Maintenance Workspace



Service Application allows your qualified service representative to adjust and optimize system parameters dependent on usage and final application.



Manager Application allows access to the system Configuration Workspace.

Note: to use this feature, it needs to be activated by your qualified service representative. It is then accessed at any time via a password (provided by the service representative). Typically, the password is kept by a user nominated Manager. The password can be changed at any time.



Consumables Application shows consumable status and allows consumable software wizards to be launched.

Note: You can visualize the same consumable status information in the **Glance Workspace > Consumables Application** but you cannot perform consumable replacement from there.

Consumables found in this application: Progard[®] pack(s) and if applicable, Vent filter and distribution loop filter.



Sanitization and Cleaning Application allows cleaning software wizards to be launched.

The RO membrane cleaning wizard will guide you through the different steps, indicate the cleaning time and what is required in order to perform RO membrane cleaning. During the first steps, until reagents are introduced into the system, you will be prompted by the wizard, to validate or cancel the RO membrane cleaning.

When opening this application, two dates related to cleanings are shown:

- The "Performed" date indicates the last time a cleaning was done,
- The "Due date" indicates when the next cleaning is due.



Note: The cleaning timers can be adjusted in collaboration with your Qualified Service representative.

The RO membrane CL_2 cleaning takes around 30 min and the RO membrane pH cleaning takes around 75 min (see Maintenance Chapter for more information).

Tank emptying: Your tank will be emptied until it reaches the 0% level of water.

TOC cleaning: TOC cleaning process will be launched, the time duration is 1 hour.

Configuration Workspace

The **Configuration Workspace** contains applications allowing modification of some system parameters

Configuration Workspace



Note: The information shown in the configuration applications can be observed, changed and saved. Configuration applications can be accessed when manager access has been activated by a qualified service representative. It is accessed by the manager using a password. In the **Workspace Glance** applications similar information can be seen, but not changed.



Information Application allows modification of system information.

Product information: Parameters from the manufacturer that uniquely identify the system cannot be changed but personalization of your system is possible by giving it a name, and a location for example.

Company name: Your Company name, address and contact information can be found and changed in this application.

Key contact (Add/Del/Mod): You can add, delete or modify your contacts phone number and email from this application. These key contacts can be viewed in the **Glance Workspace > Information Application** by any system user.



Processes Application

Makeup process: The Tank refill setpoint can be modified, the RO recovery setpoint can be decreased. The Tap feed pressure max. is given as information only.

ASM scheduling: If the option has been activated, ON time and cycle period can be adjusted based on the system daily usage. It is not recommended to disable the ASM UV lamp cycle.

Distribution process: Auto lab closed can be activated here. The auto lab closed must be activated to refresh the tank water if the system was not in tank filling for more than 3 days and the tank level is above the tank refill setpoint.

The "Auto recirculation duration" can be adjusted from 15 to 60 min per hour, and the Post recirculation from 30 to 60min per hour.

The "Tank flush setpoint" can be adjusted to the following range of value: 0-75%.

Options: The water sensor (or leak detector) can be activated here.

Measure outputs: Two measures can be selected from the list of parameters, these parameters can be monitored by external devices.

Units: Set your desired units here.



Alarm and Alert Application

Alarm thresholds: Thresholds corresponding to the main system parameters can be adjusted here. If a threshold has been crossed, the system will raise the corresponding alert or alarm.

Example: RO Low TDS setting can be activated in some particular tap feed water cases. If the tap feed water quality is too good, the RO membrane rejection performance maybe poor and based on the system RO recovery, the system could raise alarms related to the RO rejection, or tap feed conductivity. This parameter can be selected to suppress these alarms.

External signals: This option must be installed and activated by your service representative. Unlimited alerts or alarms can be selected from the list, these alarms can be reported to an external device.



System Settings Application

LCD: Adjust the brightness of the display.

Language: Your system language has been set by your Qualified Service representative. However, the language can be changed using this application.



Note: Pay attention to the fact that with an unfamiliar language you may encounter some difficulties to return to your original language.

Sound: Sounds can be configured here.

Alarm and Alert sounds can be selected independently. You can also decide to activate Keyboard clicks.

Date & Time: The system date and time can be changed here.



Note: Seasonal changes to the time (i.e. Daylight Savings) are automatic.

Network configuration: The LAN communication can be configured here (fixed ip address or DHCP can be selected)



Export/Import Application


The system configuration file can be exported from this application. It is highly recommended to export the system configuration file and keep it as a backup in order to keep your complete settings.


The configuration file can be imported into the system if it corresponds to the same system from which it has been exported previously.



Distribution Schedule Application allows scheduling of recirculation of the distribution process.

Behaviours: 3 distribution modes - Auto-recirculation, Recirculation and Tank Flush - can be programmed for each day of the week.

Auto-recirculation is set by default. Start and end times of the distribution modes can be defined and added to the distribution schedule .

Copy schedule:  the schedule defined for one day can be applied to other selected days of the week.

How to View Hidden Applications

On a workspace screen, some of the available applications are located on a second screen.

Procedure


If the workspace window contains more than the maximum application 9 available, you must change the page to access to the other available applications.

- On the left and right side of the application windows, two arrows are present. Click on these arrows to navigate through different application windows.



Clicking on the right arrow on Glance workspace will show the additional applications.



 **Note:** When entering in a workspace screen that contains more than one page of applications, the first page will always be displayed.

Alarms and Alerts - Definition

Notifications are displayed when the system has an Alarm or an Alert.

Type	Description	Examples
Alert	An alert is displayed in yellow meaning that maintenance is required or a non-critical event occurs. The system continues to operate.	Replace Progard® in xx days.
Alarm	An alarm is displayed in red, meaning that a problem was detected by the system. The system continues to operate.	Elix Resistivity < Set point.
Alarm Stop	An alarm is displayed in red, meaning that a critical problem was detected by the system. The corresponding system process is stopped : the system will not produce or distribute water until the problem is fixed.	Low feed water pressure.

How to Acknowledge Alarms and Alerts

Some Alarms stop the corresponding process in order to protect the system, acknowledgment of these Alarms resume the process if the cause of the alarm has been corrected. Alarm messages cannot be cleared unless the cause of the alarm has been fixed. Acknowledging Alerts clears the message for 24 hours.


Procedure

To acknowledge an alarm or alert:


1. Tap on the Alarms and alerts Notification button on the **Overview** screen.
Alarm blinking




The **Alarms and Alerts** screen is opened:

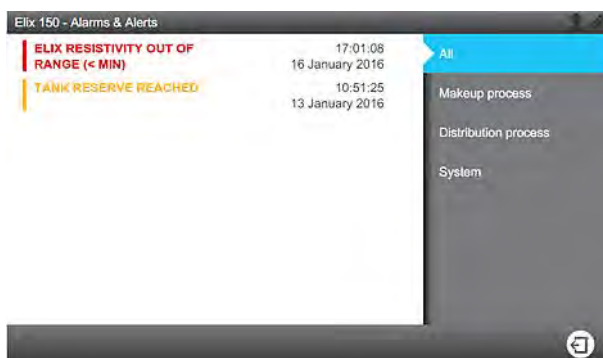

2. By clicking on a single event message a screen opens describing the alarm or alert and if required information on fixing it.

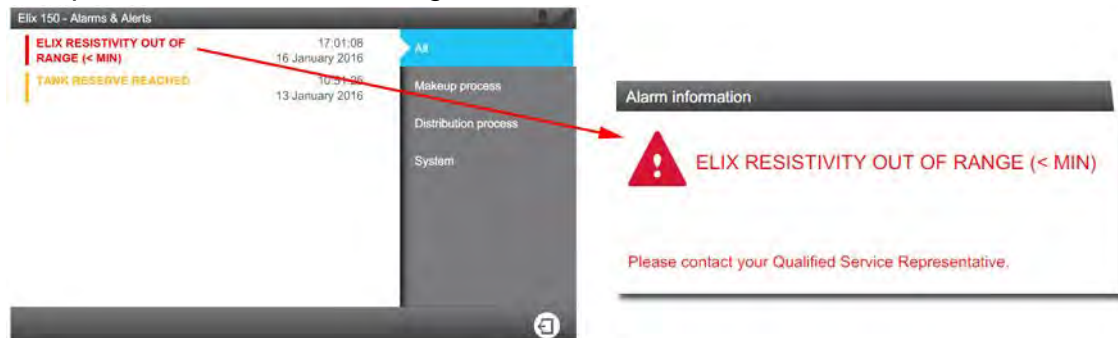

- If the alarm requires a cleaning or a consumable replacement, a software wizard opens and guides you through the process.

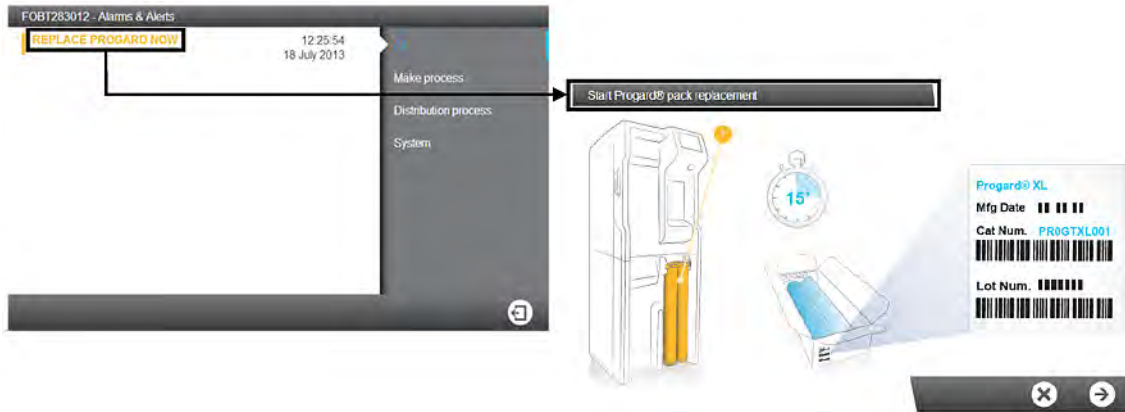

Results






- Once Alarms are acknowledged and the cause of the alarm has been fixed the process will return to Ready.

Screensaver

When the system is idle, a screensaver appears on the display.

The system continues to operate and the **screensaver** displays the main system information.

Screensaver visual states:

STANDBY		Makeup process standby and/or Distribution process are in STANDBY .
READY		Makeup process and Distribution process are in READY .
READY + ALERT		Makeup process and Distribution process are in READY with one or more alert(s) .
ALERT		One or more alert(s) are active. Maintenance to be performed .
ALARM + ALARM STOP		One or more Alarm(s) are active.

Tapping on the screen when the **screensaver** is active will wake the system and open the **Overview** screen.

Maintenance

Alerts are triggered when a Consumable Replacement, Cleaning or Sanitization is required.

How to Use Maintenance Wizard

Consumable Replacements, Cleanings and Sanitizations are undertaken using software wizards. The wizard for the action can be launched directly from the alert.

Procedure



Note: The maintenance wizards can be launched in different ways.

The maintenance wizards can be launched using the **Consumables** application:



or the **Sanitization & Cleaning** application:



See Applications at the **Maintenance Workspace** chapter:




1. Tap on the Alerts Notification button on the **Overview** screen to open the Alarms and Alerts screen.



2. By clicking on a Consumable Replacement, Cleaning or Sanitization message, a software wizard is opened to guide you through the procedure.

Start Progard® pack replacement

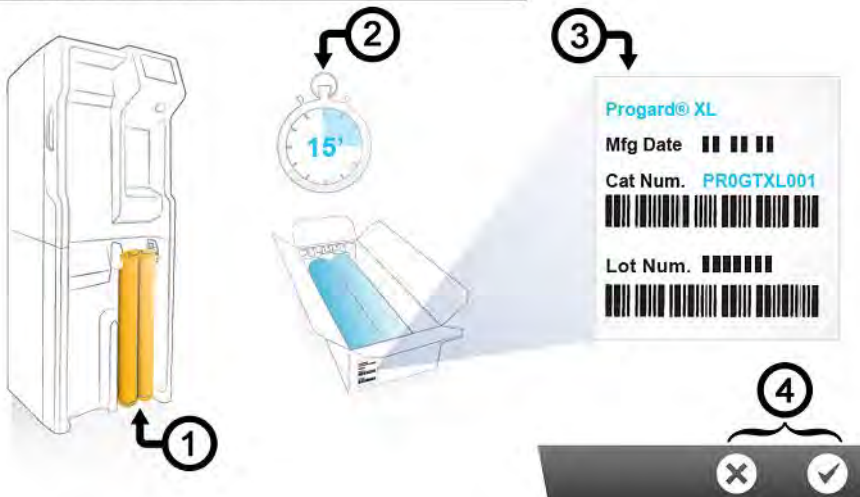

The software wizard shows:

1. The location
2. The duration (Estimated time)
3. The references needed

3. You can decide to launch or cancel the software wizard (4).

Results

Once the last step is completed the software wizard closes and the corresponding alert related to the maintenance action will disappear from the list of Alarms and Alerts.

RO Cleaning tools

What do you need to perform an RO membrane Cl₂ cleaning ?

Two solutions can be chosen to perform Cl₂ cleaning.

1. Chlorine tablets & RO regeneration tools:
 - For a system with one Progard[®] Pack: One RO Regeneration Tool
 - For a system with two Progard[®] Packs: One RO Regeneration Tool and one RO Regeneration Plug/Tool
 - Two Chlorine tablets (RO Protect C)
2. Alternatively, if you can use Progard Auto-clean Packs that already contains Cl₂ cleaning agent.

What do you need to perform an RO membrane pH cleaning ?

- For a system with one Progard[®] Pack: One Regeneration Tool
- For a system with two Progard[®] Packs: One Regeneration Tool and one Regeneration Plug/Tool
- 1 to 3 Acid or Base pH reagent pouches depending on system type (1 for HX 7040 / HR 7060, 2 for HX 7080 / HR 70120, 3 for HX 7120 / HR7170, HX 7150 / HR 7220)




Note: Your Qualified Service representative has selected, based on your tap feed water quality, the pH cleaning type that you will need to perform on your RO membranes and the pH cleaning frequency.

RO Cleaning

Table 1: RO Cleaners

RO Cleaners	Conditioning	Recommended usage	Action
Chlorine - RO Protect C (Cl ₂ tablets)	Delivered as a tablet	Use periodically when prompt by the system* every 84 days or as recommended by your Service representative.	The regular use of a chlorine tablet helps to reduce biofilm formation on the surface of the RO membrane.
RO Acid Cleaner - RO Care A	Delivered in a pouch. The empty pouch is removed after use.	Use periodically when prompt by the system* or when the RO % rejection drops more than 5% and/or when the RO permeate flowrate drops more than 10%, if due to mineral scaling.	The use of RO Acid Cleaner removes some or most mineral scale buildup on the surface of the RO membranes.
RO Base Cleaner - RO Care B	Delivered in a pouch. The empty pouch is removed after use.	Use periodically when prompt by the system* or when the RO % rejection drops more than 5% and/or when the RO permeate flowrate drops more than 10%, if due to organic fouling.	The use of RO Base Cleaner removes some or most organic material buildup on the surface of the RO membranes.

(*): RO cleaning alerts can be activated and timers adjusted based on feed water quality.

 **Important:** Do not replace RO cleaners with other chemicals. The chemical concentration and form are not designed to fit with our cleaning program, and it may damage the RO membranes & water system.

Communication Ports & Software

USB

The water system has a built-in USB port that offers the possibility to download the water system user manual or to export the system data and/or history.

USB port is located just under the main display. This is a hot pluggable port that automatically detects a USB key when a compliant device is connected:

- USB 2.0 compliant,
- Type A,
- FAT16 / FAT32 under Windows® Operating System, and
- ext3 / ext4 under Linux Operating System.

Ethernet

The water system has a built-in Ethernet port that offers the possibility to connect to a TCP/IP network.

Up to three users can connect to the system at the same time but only one single user can access applications that modify system parameters (applications found in **Configuration Workspace** as an example).

You can access to the same applications and screen views on the remote display interface and the system's main display.



Note: When three users are connected, the system informs any new connection that the maximum number of opened sessions has been reached.

Supported browsers

When connected through Ethernet protocol, the display interface can be accessed remotely using the following internet web browser versions.

The remote display of the main display is compatible with the following browser versions:

Table 2: Internet Browser Compatibility

Browser type	Recommended version
Chrome® Software	39.0.2171.71

System software

System software included in this product contains copyrighted software that is licensed under the GNU GPL.

Details of the legal notices relating the system software licence, claim, open source and components references can be found from the system MMI at **Glance / User Manual / Legal notices**.

You may obtain the complete Corresponding Source code from us for a period of three years after our last shipment of this product, by submitting a written request to your service representative.

How to Change the Network Configuration

To change your Internet or Network configuration use the System Settings Application.

About this task

! **Important:** The network configuration should only be changed from the system display. Do not attempt to change the settings remotely using the network.


Procedure


1.

From the Configuration Workspace , select the System


Settings Application .

2. Select **Network configuration**.

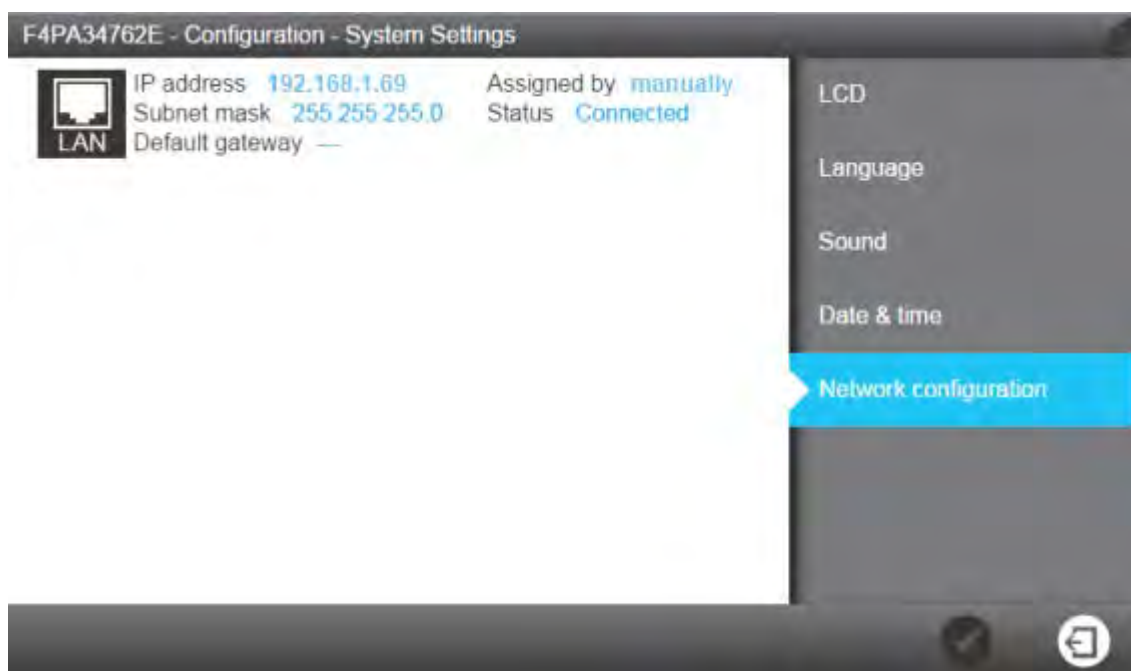

 **Note:** This screen shows the default factory settings from the system LAN configuration.

3. Click on the LAN area of the screen to open the complete LAN configuration settings.

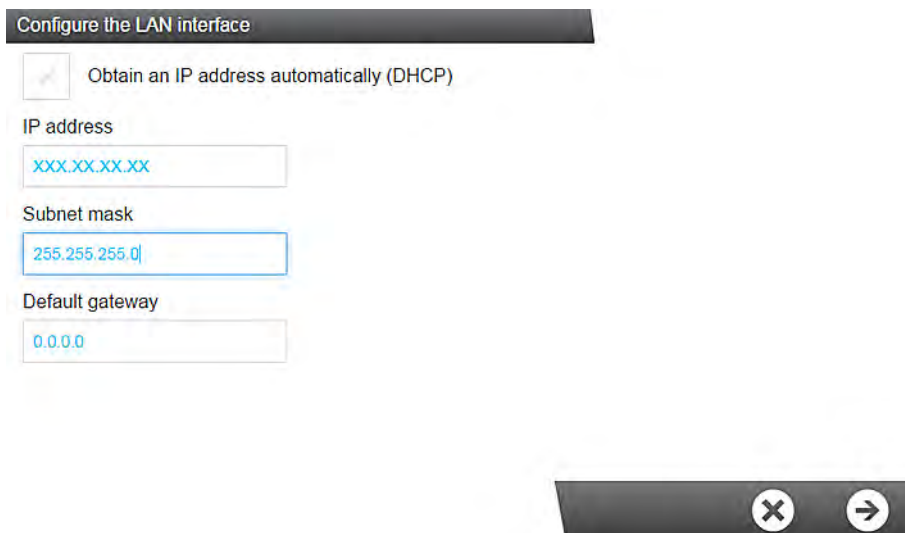


4. Enter the parameters corresponding to your configuration.


If you use DHCP mode, check the DHCP box. If using static mode you must complete the three fields:

- IP address
- Subnet mask, and
- Default gateway

using ip V4 format.

- 5.

Validate your configuration by pressing the tick  icon.

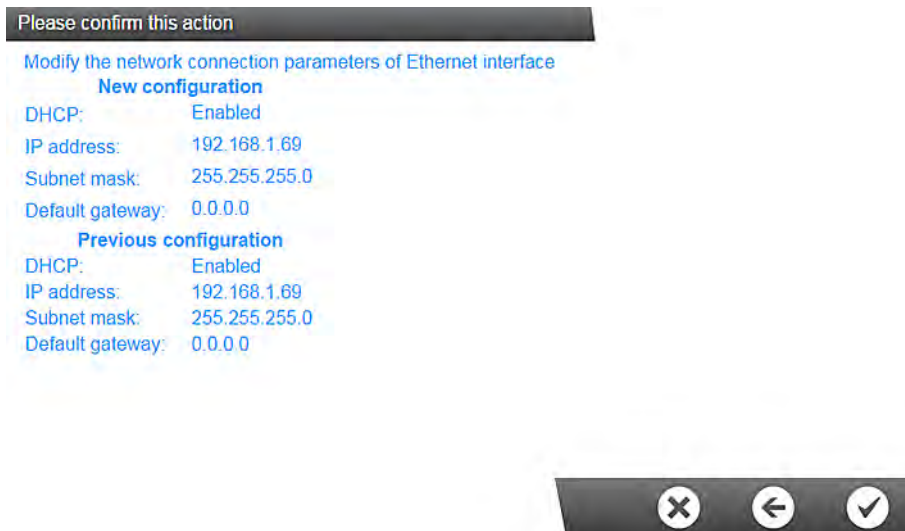


6. The Network configuration is changed and the display returns to the System Settings Application.

How to Power On and Power Off the Water System

Procedure

The water system has been designed to remain powered, this ensures the water quality in the system is maintained.

Water system backpanel view: (1) Power switch; (2) Power Outlet

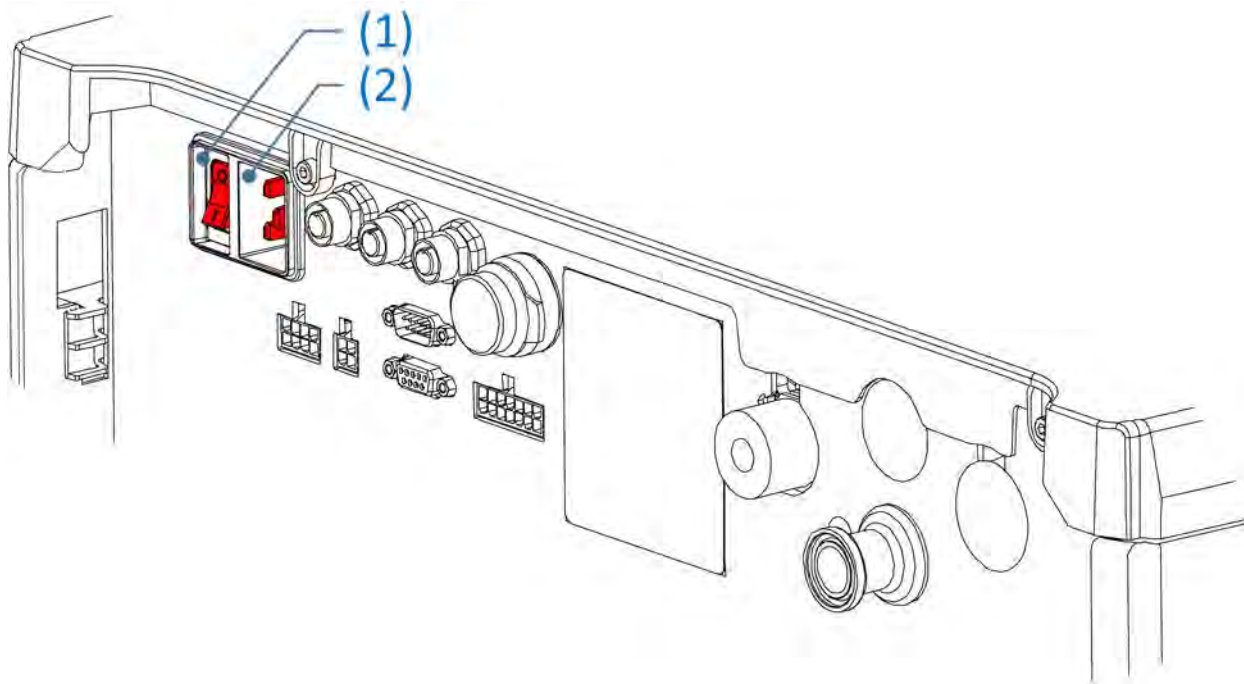

Note: Do not power off the system by removing the power cord from the **Power Outlet**.

If you need to power off the water system for any reason:

1. Put the Makeup and Distribution processes into **Standby**.
It is not recommended to shut down a running water system when operating (producing or dispensing water).
2. Use the **Power Switch** to power off the water system.

How to Use the Distribution Schedule

Procedure

The **Distribution Schedule** application from the **Configuration** workspace allow you to program the distribution loop to fit your needs.


Depending on your production or system usage, you can affect a distribution behavior to each 1/2 hour of each day of a the week. The three distribution possible behavior are: Auto-recirculation, Recirculation and Tank Flush.

The distribution schedule can be configured, within 30 minute slots, to be in Auto-recirculation, Recirculation or Tank Flush.

1. From the **Configuration** workspace, open the **Distribution Schedule** application:



2. The **Distribution Schedule** application opens (*the bullet points in red will guide you throught the following steps*):

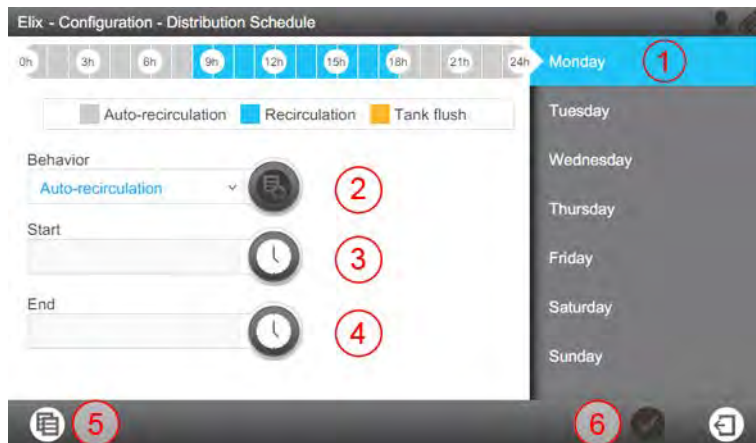

- a) Select one day of the week from the right side screen area (1).
- b) Choose the desired Distribution Schedule "Behavior" using a drop down list
- c) Program a start time to this behavior from field (3).
- d) Then program a end time (4) to this behavior.
- e) Each behavior has a start and a stop time which you can choose. After choosing these items, press the '+' button (2) to confirm your choices. You can see the result of this action at the top of the screen (0h .. 24h).
- f) When at least one slot of time has been completed, the save function (6) becomes white and allow you to save the scheduled slot of time.

- g) Once a day has been completed, a "profile" day can be copied to another day of the week using the copy function (5).

Copy

From
Monday

To

<input checked="" type="checkbox"/> Monday	<input checked="" type="checkbox"/> Friday
<input checked="" type="checkbox"/> Tuesday	<input checked="" type="checkbox"/> Saturday
<input checked="" type="checkbox"/> Wednesday	<input checked="" type="checkbox"/> Sunday
<input checked="" type="checkbox"/> Thursday	

X →


How to Flush the Tank

Procedure

This procedure can be used only if the system is managing the distribution loop.

Prerequisite: Put the Distribution processes into Standby (makeup process will be automatically stopped by the application)

1. From the **Maintenance** workspace, launch the **Sanitization & Cleaning** application and select Tank emptying (1).

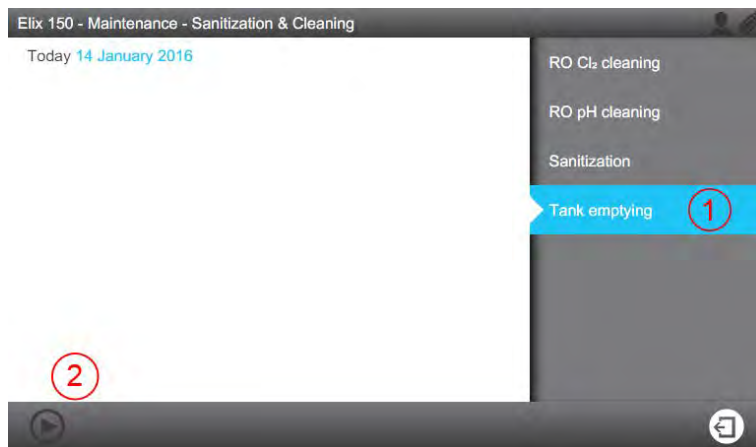

2. Launch the tank emptying function using the bottom left button (2).



Note: If the tank emptying button is disabled (grey), check if makeup and distribution processes are in Standby.

Ordering Information

Consumable Catalogue Numbers

Packs and Filters

Label	Catalogue Number	Description
Progard [®] XL-S-CL	PR0GTXLCS1	Autoclean Qty 1
Progard [®] XL-S-CL	PR0GTXLCS2	Autoclean Qty 2
Progard [®] XL	PR0GTXL001	Qty 1 (US, Canada and Mexico only)
Progard [®] XL-S	PR0GTXL0S1	Silver Impregnated Carbon - Qty 1
	PR0GTXL0S2	Silver Impregnated Carbon - Qty 2
Progard [®] XL-S		

Cleaners

Label	Catalogue Number	Description
RO Protect C	ZWCL01F50	CL ₂ tablets - Qty. 50
RO Protect C – USA only	5874316024	CL ₂ tablets - Qty. 24
RO Protect C – Canada only	5874316024C	CL ₂ tablets - Qty. 24
RO Care A	ZWACID012	RO Acid Cleaner - Qty. 12
RO Care B	ZWBASE012	RO Base Cleaner - Qty. 12

Qty 1 refers to 1 unit per box.

Saniclean paks (Required tools for cleanings)

Label	Catalogue Number	Description
RO Regeneration Tool	ZLXLCLPAK	Required by all systems.
RO Regeneration plug/tool	ZLXLPLUGP	Required also for systems with two Progard packs.

Accessory Catalogue Numbers

Designation	Catalogue Number	Description
Water Sensor (Main)	TANKLKXL1	Water sensor to be connected on system
Water Sensor	TANKLK002	Water sensor to be connected to other water sensors (Up to 3 water sensors connected in series)
External Valve	ZLXL00ESV	External solenoid valve
Degasser Kit	ZLXLDEGK2	Internal system option
External Pretreatment Cable	ZLXLPTCAB	Communication cable
Air Gap 2 inlets	AIRGAPXL2	For tubing internal diameter 10mm
Flow switch for pretreatment unit	ZLXLPTFSW	Flow drain switch for external pretreatment unit
Alarm Report Cable	ZLXLALCAB	Reports x2 alarms and x2 4-20 mA measures
External Pressure Regulator	ZLXL000PR	Feed water pressure regulator (0-25 bars)
UF Pretreatment Unit	ZUFPREUN0	Ultrafiltration filter Unit 3/4"
UF Installation kit	ZUFPREUN8	Installation kit for the UF Pretreatment Unit 3/4"
Large ELIX/RIOS™ Connection box	ZLXLCB001	Connection box
Large Elix/Rios™ Sub-D Chaining cable	ZLXLSDC01	Chaining cable (Applicable to HX and HR systems)
Resistivity kit high throughput	ZKITRES00	Resistivity and temperature measurement on distribution.
Resistivity kit & boost high throughput	ZKITRES01	Resistivity and temperature measurement on distribution including resistivity boost.
Resistivity & TOC kit high throughput	ZKITRESTOC	Resistivity, temperature and TOC measurement on distribution.
Resistivity boost & TOC high throughput	ZKITRES01T0C	Booster media and TOC measurement on distribution.






















System Catalogue Numbers


ZLXL72040	Milli-Q [®] HX 7040 Water Purification System (HC) 100V 50/60 Hz
ZLXL62040	Milli-Q [®] HX 7040 Water Purification System (HC) 120V 60 Hz
ZLXL52040	Milli-Q [®] HX 7040 Water Purification System (HC) 230V 50/60 Hz
ZLXL71040	Milli-Q [®] HX 7040 Water Purification System (LC) 100V 50/60 Hz
ZLXL61040	Milli-Q [®] HX 7040 Water Purification System (LC) 120V 60 Hz
ZLXL51040	Milli-Q [®] HX 7040 Water Purification System (LC) 230V 50/60 Hz
ZLXL72080	Milli-Q [®] HX 7080 Water Purification System (HC) 100V 50/60 Hz
ZLXL62080	Milli-Q [®] HX 7080 Water Purification System (HC) 120V 60 Hz
ZLXL52080	Milli-Q [®] HX 7080 Water Purification System (HC) 230V 50/60 Hz
ZLXL71080	Milli-Q [®] HX 7080 Water Purification System (LC) 100V 50/60 Hz
ZLXL61080	Milli-Q [®] HX 7080 Water Purification System (LC) 120V 60 Hz
ZLXL51080	Milli-Q [®] HX 7080 Water Purification System (LC) 230V 50/60 Hz
ZLXL72120	Milli-Q [®] HX 7120 Water Purification System 100V 50/60 Hz
ZLXL62120	Milli-Q [®] HX 7120 Water Purification System 120V 60 Hz
ZLXL52120	Milli-Q [®] HX 7120 Water Purification System 230V 50/60 Hz
ZLXL72150	Milli-Q [®] HX 7150 Water Purification System 100V 50/60 Hz
ZLXL62150	Milli-Q [®] HX 7150 Water Purification System 120V 60 Hz
ZLXL52150	Milli-Q [®] HX 7150 Water Purification System 230V 50/60 Hz

ZR0L72060	Milli-Q [®] HR 7060 Water Purification System (HC) 100V 50/60 Hz
ZR0L62060	Milli-Q [®] HR 7060 Water Purification System (HC) 120V 60 Hz
ZR0L52060	Milli-Q [®] HR 7060 Water Purification System (HC) 230V 50/60 Hz
ZR0L71060	Milli-Q [®] HR 7060 Water Purification System (LC) 100V 50/60 Hz
ZR0L61060	Milli-Q [®] HR 7060 Water Purification System (LC) 120V 60 Hz
ZR0L51060	Milli-Q [®] HR 7060 Water Purification System (LC) 230V 50/60 Hz
ZR0L72120	Milli-Q [®] HR 7120 Water Purification System (HC) 100V 50/60 Hz
ZR0L62120	Milli-Q [®] HR 7120 Water Purification System (HC) 120V 60 Hz
ZR0L52120	Milli-Q [®] HR 7120 Water Purification System (HC) 230V 50/60 Hz
ZR0L71120	Milli-Q [®] HR 7120 Water Purification System (LC) 100V 50/60 Hz
ZR0L61120	Milli-Q [®] HR 7120 Water Purification System (LC) 120V 60 Hz
ZR0L51120	Milli-Q [®] HR 7120 Water Purification System (LC) 230V 50/60 Hz
ZR0L72170	Milli-Q [®] HR 7170 Water Purification System 100V 50/60 Hz
ZR0L62170	Milli-Q [®] HR 7170 Water Purification System 120V 60 Hz
ZR0L52170	Milli-Q [®] HR 7170 Water Purification System 230V 50/60 Hz
ZR0L72220	Milli-Q [®] HR 7220 Water Purification System 100V 50/60 Hz
ZR0L62220	Milli-Q [®] HR 7220 Water Purification System 120V 60 Hz
ZR0L52220	Milli-Q [®] HR 7220 Water Purification System 230V 50/60 Hz




Appendix

Display Icon Description

Icon	Function
	Exits the current application or wizard.
	Navigates back to the previous screen.
	Navigates forward to the next screen.
	Cancels an action.
	Confirms an action.
	Adds a new item to a list.
	Removes the selected item(s) from a list.
	Modifies the selected item(s) of a list.
	Opens the Glance Workspace.
	Opens the Maintenance Workspace.
	Opens the Configuration Workspace.
	Shows water quality measurements.
	Shows hydraulic parameters.
	Shows electrical parameters.
	Starts a software wizard.
	Opens a consumable override wizard.
	Filters the data.
	Exports the data.
	Automatically fills a field with the system date.
	Time setting
	Copy schedule

Icon	Function
	Add selection to schedule


The system buttons used are virtual icons on the display and their **status** is determined by colour.

	Disabled.
	Enabled.
	Pressed or selected.



Peripherals and communication indicators:

On each of the MMI screens, on the top right, there are two icons to indicate the connection status, via Ethernet or the front side USB port.

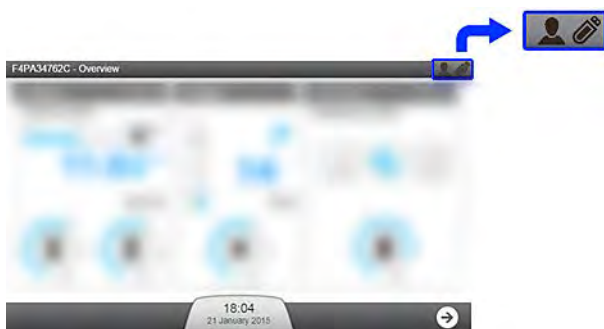

As an example for the Overview screen:



USB port (front side):

	No USB device detected on the system.
	A USB device has been inserted and detected by the system.

Ethernet connection status:

	No effective network connection has been established.
	An Ethernet connection is active on the system. At this time, up to 3 different IP address can be connected to view the general operation of the System.

 10.150.65.3

Via the Ethernet connection, someone is remotely using the **Maintenance** and/or **Configuration** applications on the System. The IP address is shown. At this time, no one else can access the **Maintenance** and/or **Configuration** applications. To gain access at this time via the Ethernet connection, ask the user having the remote IP address to leave the **Maintenance** and/or **Configuration** applications.

 Local

A user is using the MMI directly and is in the **Maintenance** and/or **Configuration** applications. At this time, no one else can access the **Maintenance** and/or **Configuration** applications. To gain access at this time, via the Ethernet connection, ask the user (in front of the system) to leave the **Maintenance** and/or **Configuration** applications.

System Operating modes

This chapter describes the various process states of Makeup and Distribution processes.

When not in standby, the different modes possible within the Makeup and Distribution processes are described below.

Table 3: Makeup Ready Modes

Makeup Mode	Use
INITIALIZATION	To check and reset components.
READY	Pauses the Makeup process when the tank is full.
FLUSH	To periodically sweep away contaminants on the feedwater surface of the RO Membrane.
RINSE	To prevent poor quality water reaching the Elix [®] Module before TANK FILLING.
TANK FILLING	To fill the Tank.
AUTOTEST	Checks internal components.
PROGARD FLUSH	To flush new Progard [®] Pack(s).
RO RINSE	To rinse new RO Membrane(s).
CL ₂ CLEANING	To clean RO Membrane(s).
pH CLEANING	To clean RO Membrane(s)
REAGENT REMOVAL	To remove cleaning agents after CL ₂ CLEANING or pH CLEANING.
ALARM STOP	Stops the Makeup process in the event of an Alarm Stop signal.
BACKWASH FILTER REGENERATION	Pauses the system whilst a backwash filter regenerates.

Table 4: Distribution Ready Modes

Distribution Mode	Use
READY	To get the distribution process operating in one of the different modes: Recirculation, auto-recirculation or schedule.
RECIRCULATION	To maintain water quality.
TANK FLUSH	To maintain water quality when tank is full.
ALARM STOP	To stop the Distribution process in the event of an alarm stop signal.
TANK EMPTYING	To empty the tank.



Note: When the system is powered on, it returns to the mode operating when the system was powered off. For example, if the system was in Makeup **TANK FILLING** then it returns in Makeup **TANK FILLING**; and in Distribution **RECIRCULATION** then it restarts in Distribution **RECIRCULATION**.

Emergency Procedure in Multi-system configuration


When your primary system is not operational, the first auxiliary system can be converted to a primary to continue production of purified water.

Procedure

! **Important:** This emergency procedure is as temporary way to ensure continuity of production. During the emergency procedure, the consumables and actuators lifetime on the Distribution part are not managed by a system. Call your Service representative to get your primary system up and running.

1. Power off the primary system using the main switch on the system backpanel.
2. From the primary system backpanel, disconnect the three external connectors (C1, C2 Out and CH In) . The cable CH Out can be left in place.

Note: If the distribution is not duplex, C2 Out is not present.

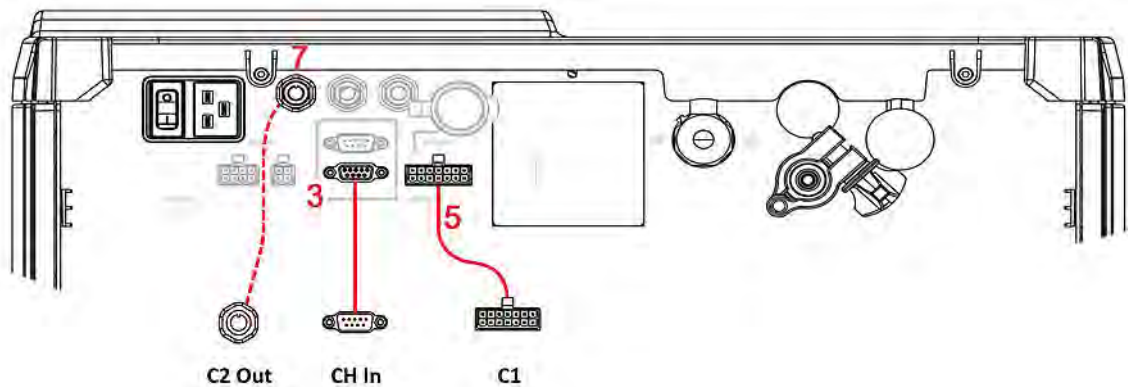

3. Connect the three external connectors to the backpanel of the first auxiliary system I. The first auxiliary will become the temporary primary system.
4. On the first auxiliary system, from the Overview screen on the MMI, press the three horizontal bars on the distribution top right panel to activate the distribution.

Before



After

